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Small size and strong power

Powered by **FITTOP**

WSHMINI

Brugermanual: 3-15
Brukermanual: 17-29
Bruksanvisning: 31-43
Käyttöohje: 45-57
User Guide: 59-71

INDHOLD

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Important information

STATEMENT

To reduce the risk of malfunction, electric shock, injury, fire, death, and equipment damage, please read all instructions and warnings in the product user manual before using.

Prohibition

1. In case of muscle injury please consult your doctor before using product. Minors are advised to use product under the adult supervision.
2. Do not use on wet and dirty skin. It can be used over the clothes, just gently press and move on the skin for about 60 seconds on each part
3. Do not use over the neck, bones and body parts with less muscle coverage. It's only applicable on the body soft tissue without causing pain or discomfort.
4. Do not use when the body has bruise, pain and discomfort.
5. Do not put any objects in the vent holes of the product.
6. Do not immerse the product in water and please keep water out of the product vents.
7. Do not put the product into the fire or close to the fire source.
8. Do not drop this product intentionally.

9. Do not disassemble this product by any means.
10. Do not close to the axis and back of the massage head by fingers, hair or other parts of the body, as it may be clipped.

Do not touch this part when starting the product without massage head installed.

Do not touch this part when starting the product with massage head installed.

Warnings and reminders


1. Keep away from liquids or heat sources and keep the vents dust free.
2. After the charging is completed, pull out the charger.
3. The product cannot be used during charging.
4. Please use this product according to the instructions.

5. Carefully check whether the product can be used normally according to the instructions in this book before each use.

Medical note:

The product cannot be used in any of the following cases without the prior consent of the doctor.

Pregnancy, diabetic complications (such as neuropathy or retinal damage), wearing a pacemaker, recent surgery, epilepsy or migraine, disc herniation, spondylitis, vertebral detachment, or spinal joint stiffness, recent joint replacement surgery or IUD, metal needle or steel plate, and any other serious problems related to your health. These contraindications do not mean that you cannot use this product. We recommend that you consult a doctor before use.

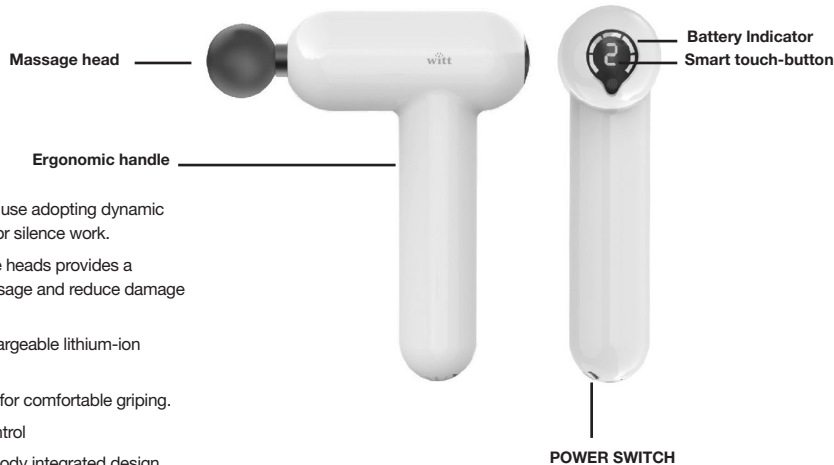- Mute technology use adopting dynamic force balancing for silence work.
- Silicone massage heads provides a comfortable massage and reduce damage to bones.
- High quality rechargeable lithium-ion battery.
- Compact handle for comfortable gripping.
- Touch screen control
- Battery and the body integrated design gives ergonomic handle with no risk of losing it.

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Massage heads installation method:

While pulling/installing, slightly rotate the massage head to pull out/install. Do not pull/install directly. Choose massage head position for most comfortable massage.

Basic parameter

Product name: Witt Super Hit Mini	Rated Voltage: 5V
Product model: WSHMINI	Rated current: ca. 2000 mA
Product weight: ca. 0,6 kg	Massage heads: 3 different massage heads
Power supply: rechargeable lithium battery	Mode: 5 speeds optional, manual speed control

Charging

1. Charge the battery before initial use. (If the product is sufficiently charged, there's no need for charging).
2. When charging, insert the USB cable into the charger plug and connect it to the general household power outlet and connect the other end to the device charging port.
3. Battery level displays when charging: battery power <.25% (red light flashing), battery power <50% (yellow light flashing), battery power <75% (cyan light flashing), battery power >75% (green light flashing), battery power 100%(while charging finishes, green light is always on).
4. The device cannot be started up while charging.

Instructions

1. Turn the power switch to the ON position, the machine stands by.
2. Speed mode: Short press the button  to turn on, 1st speed by default; short press the button  again to enter 2nd speed; short press the button  again to enter 3rd speed; short press the button  again to enter 4th speed; short press the button  again to enter 5th speed; then short press the button  again to enter standby status.
3. Manual speed control mode: In any mode, long press the button  for 0.5 seconds to enter the manual speed control mode, the speed increases with the time the button  is pressed (the same goes for the power indicator) until it reaches the maximum speed. Short press the button  to enter status. When not using this device, please make sure the power switch is set to the OFF position, otherwise, the battery will run out.

Maintenance, Cleaning, Storage, Transport

1. Clean the device with a slightly wet towel and dry with a soft cloth.
2. Place the product in the package provided while storing or travelling.
3. Stored in a dry, well ventilated place with non-corrosive gas and no less than 250 mm height from the ground.
The number of stacking layers does not exceed 5 layers.
4. If not in use for a long time, store in the package after cleaning and charge every six months.
5. This product is suitable for car, train, ship and aircraft transport.
6. This product must be fixed during the transit to prevent vibration, collision, roll and fall and place it upward, lightly to somewhere rainproof.

UK

Distinguish from traditional fascia therapy device

Witt Super Hit™ Mini has each massage head redesigned

Mini general massage head

Mini general masage head for better fascia stimulation, suitable for multiple muscle groups! General mode suitable for all fascia parts.

Trigger point massage head

Mini trigger points massage head for deeper penetration into the “trigger points” and resolving these “adhesive pain points” one by one!

Dedicated Trigger points release mode for whole body trigger points relieving.

Achilles massage head

Mini achilles massage head fit adequately for ankles and achilles tendon, widely affected by exercisers and practitioners. Well performed in bar muscles massage like hind neck and tendon.

GUARANTEE

A two year guarantee on manufacture and material defects is granted on your new appliance valid from the documented purchase date. The guarantee includes materials, salary and transport. When contacting service please state the appliance model and serial number.

The guarantee does not cover:

- Defects and damages which are not caused by manufacture and material defects.
- If original spare parts have not been used.
- If the instructions in the manual have not been carried out.
- If the installation has not been carried out according to this manual.
- If unauthorized personnel has carried out installation or repair work.

Transport damages

A transport damage observed when the dealer delivers the appliance at the customer is solely a matter between the customer and the dealer. If the customer has carried out the transport alone the manufacture / supplier cannot be held responsible for transport damages. All transport damages must be reported no later than 24 hours after delivery otherwise all claims will be rejected.

Unnecessary service visits

If a service technician is called out to a defect and the customer could have corrected the defect without help e.g. by following the instructions within this manual or by changing a fuse the customer will be charged all service costs – even under guarantee.

Commercial purchase

Commercial purchases are any purchase of appliances not used in a domestic household but instead are used for commercial purposes e.g. in a restaurant, café etc. or used to rent out or used for any other purpose including several users. In case of commercial purchase no guarantee is granted on this appliance since this is only meant for domestic household use.

SERVICE

For service and spare parts in England and Ireland please call below-mentioned number and you will be directed accordingly or send us an email.

Witt UK & Ireland Ltd.

(44) 148 351 6633 / info@witt.dk

WARNING: Witt Super Hit™ Master and Witt Super Hit™ Mini are not designed or intended to be used as medical equipment or to replace medical equipment. The products are not intended to prevent, diagnose or treat any illness or condition, even when used as instructed. The products are not medical equipment and are not intended to be used to diagnose illnesses or other conditions or for healing, mitigating, treating or preventing illnesses. Always contact your doctor if you have any medical questions.

The Witt Super Hit™ Master and Witt Super Hit™ Mini massage pistols are intended to ease bodily stress and pain by manipulating the body's soft tissue. Always contact your doctor immediately if you have any medical questions.